

Dental Project di Augusto Grassi

SOLUZIONE OVERDENTURE

ANCORAGGIO A BARRA

CHIAVE DI CHIUSURA

Fissaggio con gesso per impronta dei transfert usati per l'impronta di posizione in silicone.
Questa chiave di chiusura serve per verificare la posizione corretta dei transfert.
In caso di rottura del gesso durante la prova si deve fissare nuovamente la zona rottura con del gesso riprendendo l'impronta di posizione in silicone.





PROGETTO CAD

Dental Project di Augusto Grassi




Barra Fresata


Fresatura in CAM della barra con filettatura 1,6 mm per attacchi Ot-Equator.
Nella zona terminale della barra ho inserito due Vario soft come stabilizzatori.



Barra Fresata

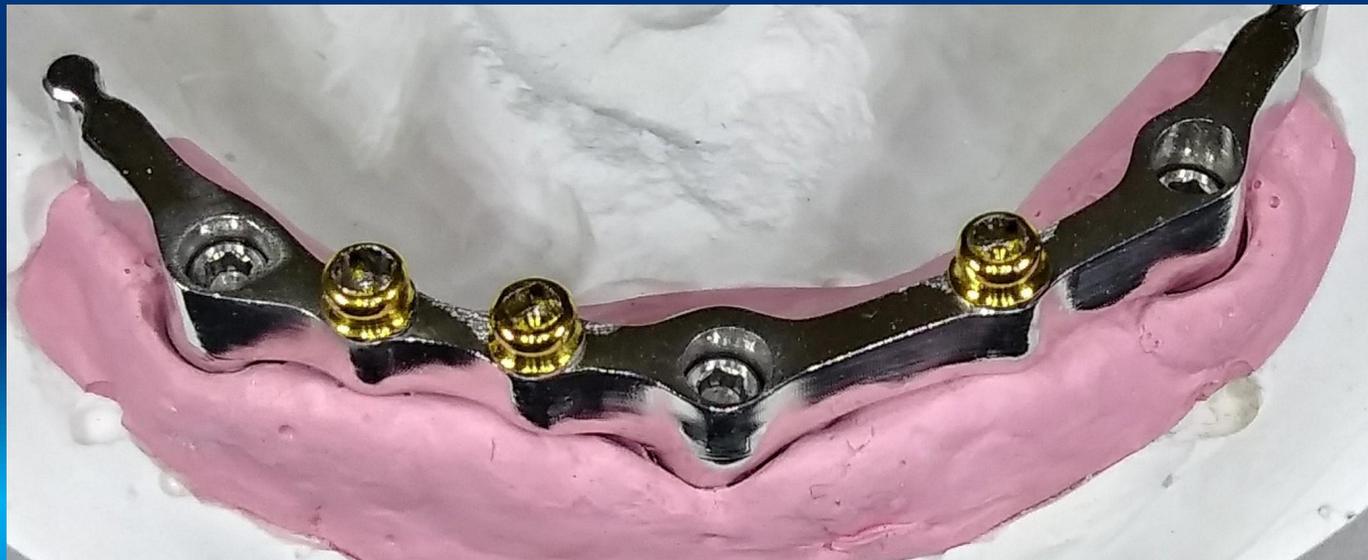


Posizionamento degli ancoraggi

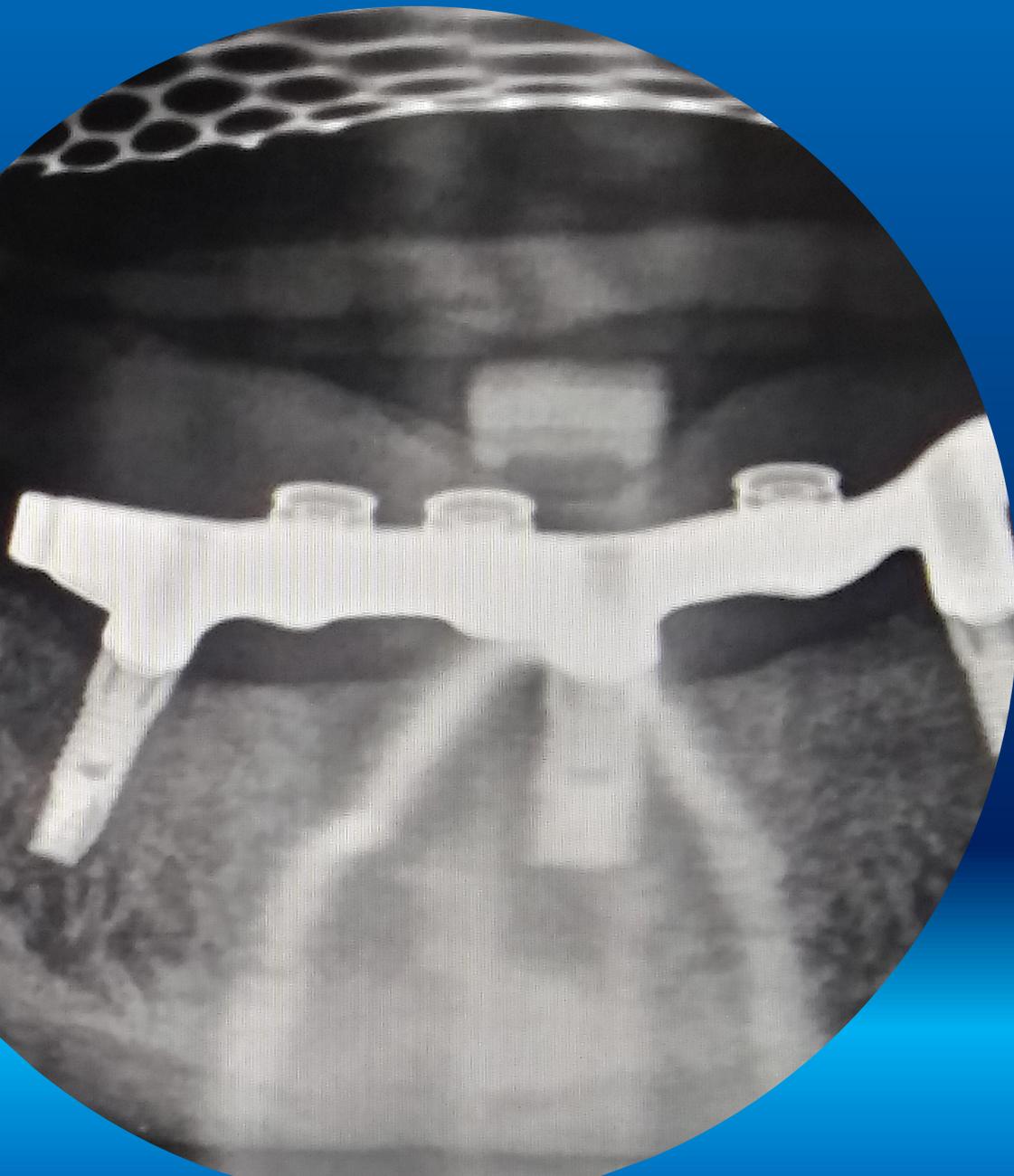

Prova Barra

RX di controllo chiusura della struttura

SOVRASTRUTTURA

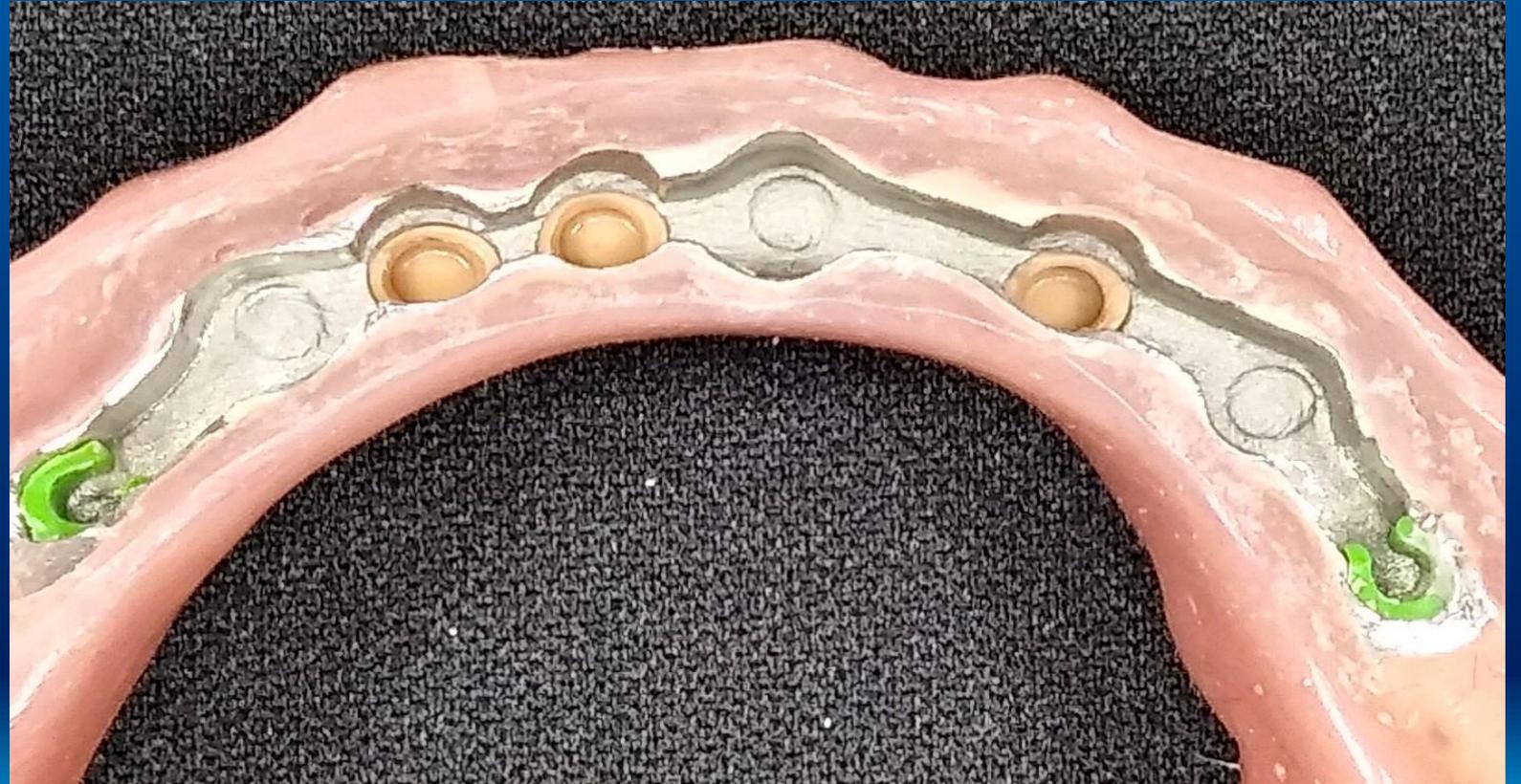



MONTAGGIO DENTI

Dental Project di Augusto Grassi



FINITURA DELLA PROTESI

Dental Project di Augusto Grassi

A close-up photograph of a row of upper teeth. A white, glossy dental material is applied to the lingual (inner) surface of the teeth, covering the root and extending slightly onto the crown. The material has a smooth, reflective texture. The background is dark, and the gum tissue is visible at the bottom.

CONCLUSIONE

Dental Project di Augusto Grassi